

**CleanSpace**<sup>®</sup>  
RESPIRATORS

EVERY BREATH COUNTS

POWERED RESPIRATORS  
FOR HEALTHCARE



CE  
2056

# A MAJOR INNOVATION IN PERSONAL RESPIRATORY PROTECTION FOR HEALTHCARE –

Designed by biomedical engineers. A lightweight compact powered respirator. Minimal parts. Ergonomic with no cables, hoses or helmets.

**CleanSpace** provides fresh clean air on-demand as you breathe and is designed to deliver high levels of protection and the ease of use required by healthcare professionals.

Delivers significant benefits in high protection, superior wearer comfort and streamlined cleaning-care protocols.

## POWER UNIT:


- Light Weight (350g/1lb)
- Airflow: 230L/min eliminates heat, moisture and fogging
- Light lithium polymer battery: up to 9hrs run time and rapid charge (<2hrs)
- Smart intuitive operating system with battery and filter alerts

## DESIGN:

- Minimal parts for simpler donning/doffing, cleaning and inventory
- Intuitive system. Fast and easy training
- Customisable coloured panels to denote personally issued devices or for use in marked zones

## MASK:

- Soft medical grade silicone. Easy and reliable cleaning
- Excellent communication. Clear voice transmission and transparent face piece.
- Low profile fits under hoods and face shields.

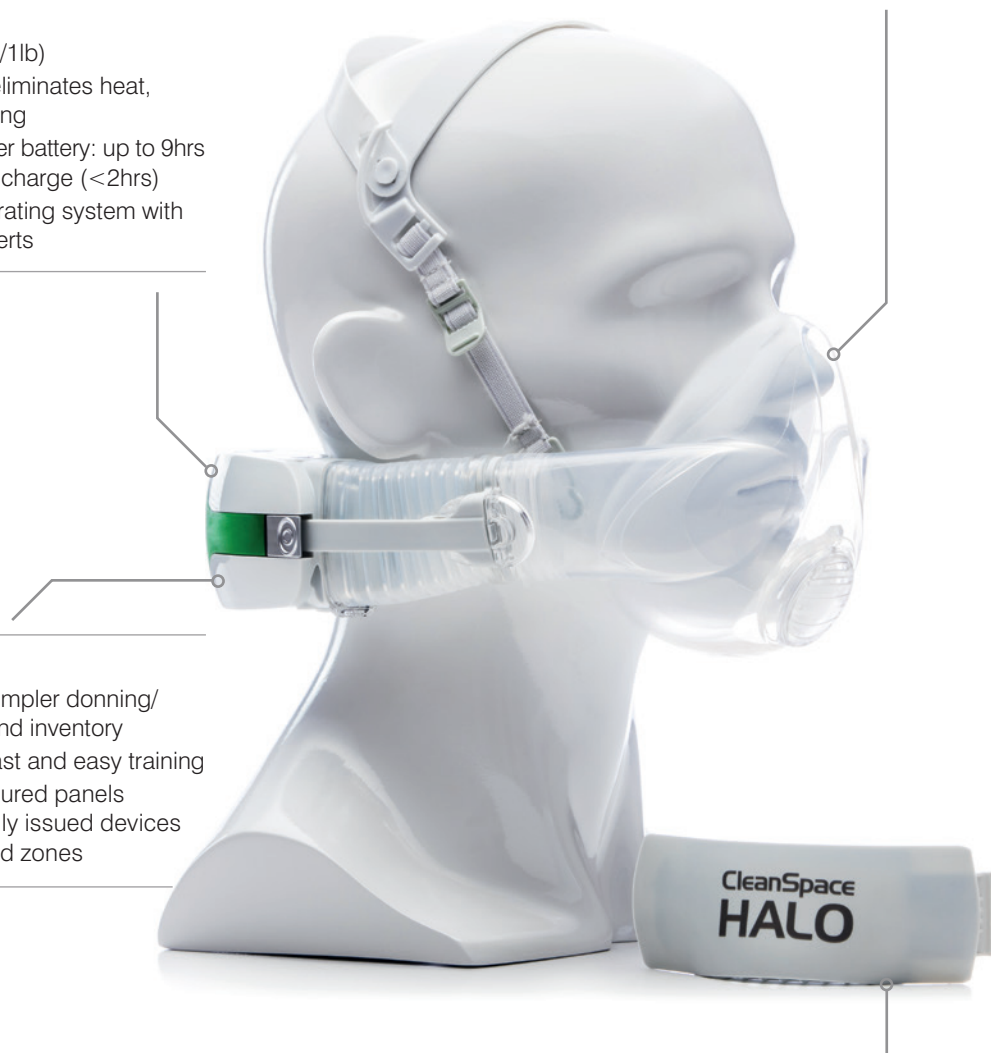

## FILTER:

- P3 protection. 99.97% filtration efficiency for 0.3micron and above