

Leading Innovation in UV-C Infection Prevention

Surfacide® provides an evidence-based, automated UV-C decontamination technology that data indicate is effective against multi-drug resistant organisms including C.diff, MRSA, VRE, CRE and Acinetobacter. The Surfacide Helios system incorporates multiple emitters that disinfect* healthcare environments in a single cycle - including the bathroom.

1 Shadows Matter

Data show UV-C is less effective on surfaces in shadowed areas. Surfacide's multiple emitters operating together significantly reduce shadows and deliver energy more effectively versus first generation UV robots.

2 Greater Coverage


Distance is critical to UV-C efficacy. As the distance increases, less energy is delivered to surfaces. Surfacide positions multiple emitters closer to high-touch surfaces, delivering more energy throughout the room in a single cycle.

3 Decreasing Time & Labor

Single emitter systems should require frequent repositioning and extended run times. Surfacide's automated system provides shorter cycle times and increased staff productivity without repositioning.

Disinfection or *disinfect is defined as the reduction of pathogenic bioburden.

surfacide[®]
UNMATCHED VISION

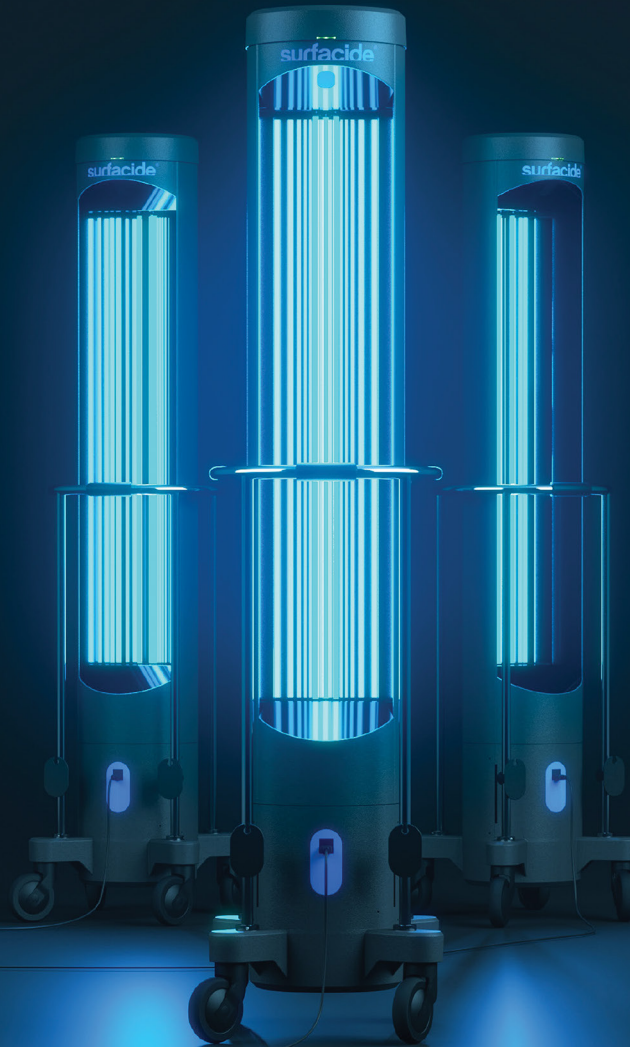

ONE
IS NOT ENOUGH

Why Surfacide?

Create a healthier environment

Join more than 350 leading hospitals around the world that have completed 1 million+ Surfacide cycles.

Unparalleled Innovation

Leading the UVC market with comprehensive advancements providing you with the power to do more.

Meaningful Design

Easy-to-use technology purposefully engineered from top to bottom for unmatched performance and mobility.

Flexibility

Emitters operate individually or in multiples according to the size of the space. Unique "Scrub" mode targets defined areas unable to be addressed by single emitters.

SURFcare

Hassle-free support including over-the-air software updates, accessory and equipment exchanges, and cloud-based portal for peace of mind.

SURFcloud

Real-time data archive provides actionable reports to drive outcomes and utilization.

Training and Support

Multidisciplinary team of experts focused on the goals necessary to achieve your success. Onsite and virtual ongoing training provided.

Cost Advantages

Greatest value in UV disinfection including 3 UV-C emitters, control tablet and multiple safety sensors all priced similar to most single emitter systems.

Multilingual

Surfacide support personnel, tablet app and manual are all available in multiple languages.

surfacide®

Helios® UV-C Disinfection System

77" emitter height
ensures disinfection
from floor to ceiling

Tower identifier

01

surfacide®
Helios® UV-C Disinfection System

Touchscreen tablet

Multilingual Surfacide App



Wireless connection
to emitters

Cloud based real-time
reporting

Joiner allows emitters
to be moved easily
in a single group

Motion sensor

CAUTION
UV DISINFECTION

surfacide

surfacide®

LED status indicator

Laser room-mapping
feature ensures

Amalgam-based
UV-C lamps

UV-C energy meter
confirms sufficient
energy output

Parabolic reflector
concentrates energy
to deliver it to where
it's needed - to all
surfaces

السرود

روداد للخدمات الطبية
ROWAD MEDICAL SERVICES

www.alrowad-ksa.com
Al sulaimaniyah Dist - Riyadh
Saudi Arabia
info@alrowad-ksa.com
mob: +966 58 274 8097