

# Homecare

## Heyer CPAP Device-A2 series



رواد للخدمات الطبية

ROWAD MEDICAL SERVICES

🌐 [www.alrowad-ksa.com](http://www.alrowad-ksa.com)


📍 Al Sulaimaniyah Dist

Riyadh, Saudi Arabia

@ [info@alrowad-ksa.com](mailto:info@alrowad-ksa.com)

☎ +966 11 400 5640

# Description

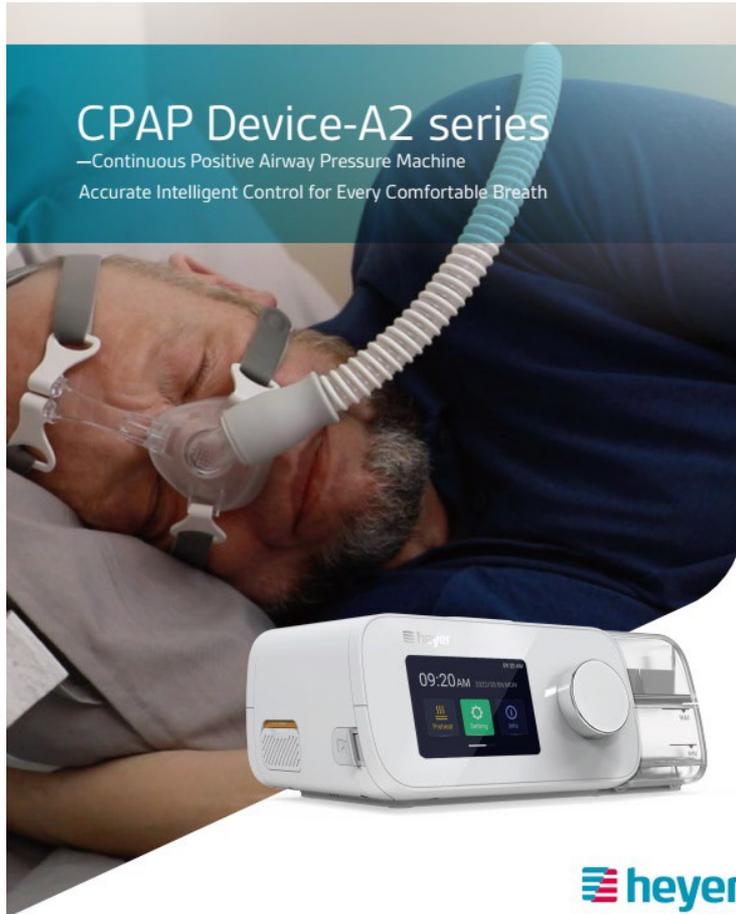

## CPAP Device-A2 series

—Continuous Positive Airway Pressure Machine

Accurate Intelligent Control for Every Comfortable Breath



## Product Specification

Product Model	20C	20A	20A Plus	20A Pro
Therapy Mode	CPAP	Auto CPAP/CPAP	Auto CPAP/CPAP	Auto CPAP/CPAP
Pressure Range (cmH2O)	4-20	4-20	4-20	4-20
Auto Ramp Time (min.)	0-60	/	0-60	0-60
Auto As-Elex Pressure Release	1-7 Level	1-7 Level	1-7 Level	1-7 Level
Identification of Central Apnea	Yes	Yes	Yes	Yes
AS-Dry Technology	Yes	Yes	Yes	Yes
Auto On/Off	Yes	Yes	Yes	Yes
Heating Tubing	Yes (Optional)			
Cellular	LTE Cat M1 & Cat NB2			
Bluetooth	Yes (Optional)			
Leaking Compensation	Yes	Yes	Yes	Yes
Titration	/	/	/	Yes
Therapy Data Storage	Yes, SD Card			
Automatic Humidification	Yes	/	Yes	Yes
Screen Size	3.5 Inches LCD Screen			
Machine Size	267 mm * 157 mm * 106 mm			
Weight	1.5 kg			
Working Noise	28±2 dB			
Power Adapter	AC Input: 100-240V, DC Output: 24V			

# Properties

**Wide-viewing  
Operation Interface**  
3.5 inches LCD Color Screen



**Updated  
Water Tank**

Big Opening and Easy for Refill  
and Cleaning  
Preheating Time up to 30 Min.



**Efficient  
Communication System**  
Cellular/Bluetooth Module/SD  
Card



**One-click Initiation  
and Cessation**  
Click Standby/Home button  
to Start/Stop Therapy




**Noise Level**  
28dB



**Multiple Choices  
of Tubing**  
5 Heating Levels  
Heating Tube with Dia.  
at 15mm  
Standard Tube with Dia.  
at 15mm or 19mm




- Adapt to a variety of masks
- Heating Tube Available

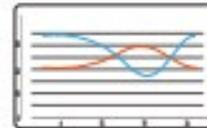


## Intelligent Temperature and Humidity Control

The CPAP system intelligently adjusts temperature and humidity levels in response to environmental variations, ensuring optimal comfort during breathing.

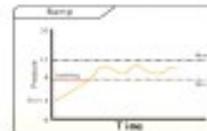

## Auto Ramp Time: Max. Up to 60min

The system intelligently detects the user's sleeping status and automatically adjusts the ramp time to align with their duration in bed, ensuring precise treatment pressure once asleep.

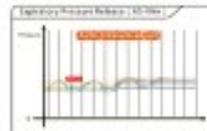

## As-Elex Pressure Release: 7 Grades x 0.5 cmH2O:

AS-Elex provides 7 grades of pressure relief during expiration, each grade with 0.5cmH2O differences. Users have more choice to adjust the pressure relief to avoid discomfort, which effectively increases patient compliance.